# **INSTRUCTIONS**

# MODEL AVX-S1-MD BIAS INSERTION UNIT

S.N.: <u>10282</u>

### <u>WARRANTY</u>

Avtech Electrosystems Ltd. warrants products of its manufacture to be free from defects in material and workmanship under conditions of normal use. If, within one year after delivery to the original owner, and after prepaid return by the original owner, this Avtech product is found to be defective, Avtech shall at its option repair or replace said defective item. This warranty does not apply to units which have been dissembled, modified or subjected to conditions exceeding the applicable specifications or ratings. This warranty is the extent of the obligation assumed by Avtech with respect to this product and no other warranty or guarantee is either expressed or implied.

## TECHNICAL SUPPORT

Phone: 613-226-5772 or 1-800-265-6681 Fax: 613-226-2802 or 1-800-561-1970


E-mail: info@avtechpulse.com World Wide Web: http://www.avtechpulse.com

# TABLE OF CONTENTS

WARRANTY2
TABLE OF CONTENTS4
FIG. 1: PULSE GENERATOR TEST ARRANGEMENT5
GENERAL OPERATING INSTRUCTIONS
AVX-S1-MD FUNCTIONAL EQUIVALENT CIRCUIT7
AVX-S SERIES DATA SHEET8
AVX-S SERIES SPECIFICATIONS

Manual Reference: /fileserver1/officefiles/instructword/avx-s/sn10282,AVX-S1-MD.doc, created July 31, 2002

### FIG. 1: PULSE GENERATOR TEST ARRANGEMENT

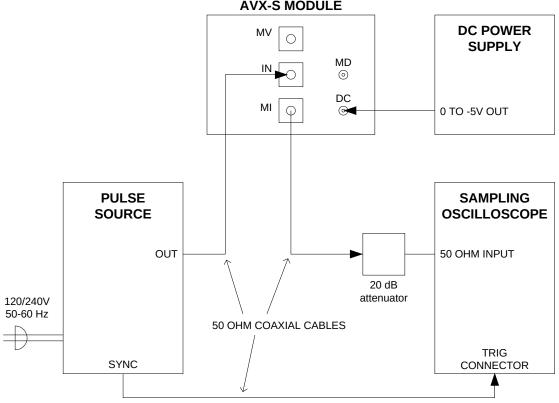

**AVX-S MODULE** 

### **GENERAL OPERATING INSTRUCTIONS**


- 1) A general description of the AVX-S1 module is given in the enclosed data sheet.
- 2) The AVX-S1 module should be connected to your pulse source via a 50-Ohm cable (supplied).
- 3) The laser diode plugs directly into the socket on the side of the AVX-S1 module.
- 4) A forward DC bias <u>may</u> be applied to the laser diode by connecting a DC potential of 0 to -5 Volts to the DC solder terminal. Note that the DC solder terminal <u>must</u> be shorted to ground if the DC bias is not used.
  - 5) The diode pulse current may be monitored by connecting the MI and MV output ports to the sampling scope via 20 dB attenuators. The output amplitude (V<sub>MI</sub> and V<sub>MV</sub>, Volts) and diode current (I<sub>D</sub>, Amp) are related as follows:

$$I_{D} = 0.2 (V_{MI} - V_{MV})$$

The laser diode voltage is given by the following:

$$V_D = 10 V_{MV}$$

#### AVX-S1-MD FUNCTIONAL EQUIVALENT CIRCUIT

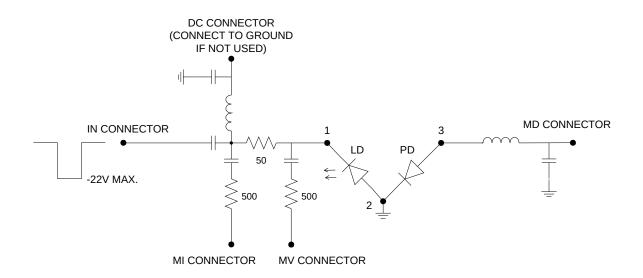

#### AVX-S1-MD FUNCTIONAL EQUIVALENT CIRCUIT (S/N 10282)



PACKAGE (S/N 10282)

## AVX-S SERIES DATA SHEET

# AVX-S SERIES SPECIFICATIONS